



Congratulations on choosing OrthoPulse® lite!



You're starting an exciting journey to a brighter, healthier smile, powered by advanced light-accelerated technology.



What is OrthoPulse® lite?

OrthoPulse® lite is an FDA-cleared medical device, featuring an exclusive, patented light therapy technology, that helps speed up your orthodontic treatment. Using **gentle near-infrared (NIR) light**, OrthoPulse® lite stimulates the bone around your teeth, making them move more efficiently. With just **5 minutes per arch per day**, you can achieve your new smile faster and more comfortably!

Getting Started with OrthoPulse® lite



STEP 1: CONNECT & POWER ON

- Connect the mouthpiece cable to OrthoPulse® lite battery.
- Short press power button to waking the device.
- LEDs indicate the battery status.



STEPS 2 & 3: INSERT & START

- 2. Place mouthpiece between teeth biting gently.
- 3. Short press power button to start the 5-min session (2 beeps and battery lights on).




STEP 4: COMPLETE SESSION

- Device beeps 3 times and status light turns off.

STEP 5: SWITCH ARCH


- Flip mouthpiece, repeat steps 2-3-4.

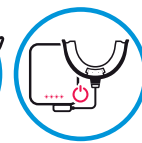



STEP 6: CLEAN & STORE

- Unplug, rinse mouthpiece, place in case to dry. Recharge battery if needed*.

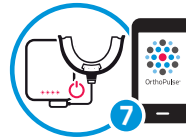

STEP 7: SYNC WITH APP (OPTIONAL)

- If paired at your doctor's office, check your SMS for the app link.
- Not paired yet? Scan the QR code to download, register, and pair your device using the serial number.
- How to sync your OrthoPulse® lite?
 - Long press power button to activate Bluetooth.
 - Light blinks while connecting. Once paired, light blinks in pairs.

How to care for your OrthoPulse® lite?

- **Keep it clean:** After each use, rinse the device under warm water and place it in the case, which allows for natural drying (no need to dry it).
- **Keep it charged and safe:** When turning on the device, if only one LED blinks (25% battery), please charge before treatment*. A full recharge takes about 6 hours. Store the device in the case.

OrthoPulse® lite Troubleshooting Tips

- **Device won't turn on?** Ensure the device is **fully charged**.
- **How to reset your OrthoPulse® lite?** Press and hold power button for 2 seconds to reset the device. You will hear 2 beeps.
- **Having trouble syncing?** Make sure **Bluetooth is on** and you're logged into the app and that it's open and updated. Sync your device data by pressing the OrthoPulse® lite power button for 2 seconds. The lights blink while connecting, once paired the lights blink in pairs.

If issues persist, contact **OrthoPulse® Support** at support@orthopulse.com or visit our website www.orthopulse.com

* Power cable must be removed to start a treatment.



M-TD-OPi-051 V1.0
OrthoPulse® lite · Model Number OPi2L-100 · UDI-Di: 91201180900032
Complies with FDA regulation and EU MDR (Regulation (EU) 2017/745)
© 2025 Biolux Technology. All rights reserved.



Scan to watch the OrthoPulse® lite Quick Start video and follow the simple steps to begin your treatment

orthopulse.com

